

DANDI:001416: Allen Institute OpenScope - Psycode (Legacy)

Region assignment: these files give no peak-channel column, and each unit lists all channels of its probe. Each unit's region is the electrode location nearest to units.depth / units.estimated_x (computed by this export). Tested on DANDI:001417, where the peak channel is given (extremum_channel_index), this method agreed on 95.6% of 3,287 units.

Scope of this export (computed by this export): the units, electrodes and intervals tables of every NWB file were read remotely from the DANDI Archive; voltage, LFP, video and other time series were not downloaded.

Source labels: **DANDI metadata** = the Dandiset's published metadata; **NWB file** = values or text stored in the NWB files; **Computed by this export** = derived by the export scripts in dandi_export/ (definitions given where used).

1. The dataset as published

Everything in this section is quoted from the Dandiset's published metadata or, for file-level fields, from sub-767925_ses-ecephys-767925-2025-02-13-11-40-19_ecephys.nwb.

Field	Value (DANDI metadata)
Name	Allen Institute OpenScope - Psycode (Legacy)
Identifier / version	DANDI:001416 / draft
Citation	Peene, Carter (2026) Allen Institute OpenScope - Psycode (Legacy) (Version draft) [Data set]. DANDI Archive. https://dandiarchive.org/dandiset/001416/draft
Contributors (roles)	Peene, Carter (ContactPerson)
License	CC-BY-4.0
Files / size / subjects	1 files, 16.80 GB, 1 subjects
Species	Mus musculus - House mouse
Approach	electrophysiological approach
Measurement technique	multi electrode extracellular electrophysiology recording technique; spike sorting technique; analytical technique; surgical technique; signal filtering technique
Variables measured	ElectrodeGroup, ProcessingModule, ElectricalSeries, Units, LFP

Description (DANDI metadata, verbatim)

Unused dandiset. For OpenScope Psycode go to dandiset 001417 <https://dandiarchive.org/dandiset/001417/draft>

File-level fields (NWB file, verbatim)

Field	Value (NWB file)
session_description	NWB file generated by AIND pipeline
institution	Allen Institute for Neural Dynamics
subject	subject_id: 767925; species: Mus musculus; genotype: Sst-IRES-Cre/wt;Ai32(RCL-ChR2(H134R)_EYFP)/wt; sex: M; age: P144D
devices	Device: Ecephys probe. Automatically generated.; ProbeA: Model: - Serial number: 20097906352; ProbeB: Model: - Serial number: 20097906092; ProbeC: Model: - Serial number: 20097909172; ProbeD: Model: - Serial number: 20097906681; ProbeE: Model: - Serial number: 20097906241; ProbeF: Model: - Serial number: 17131307182
electrode_groups	ProbeA: description = 'no description', location = 'unknown'; ProbeB: description = 'no description', location = 'unknown'; ProbeC: description = 'no description', location = 'unknown'; ProbeD: description = 'no description', location = 'unknown'; ProbeE: description = 'no description', location = 'unknown'; ProbeF: description = 'no description', location = 'unknown'

2. Contents of the files

Quantity	Value	Source
NWB files read / subjects	2 / 2	NWB file
Species	Mus musculus (2)	NWB file
Sex (files)	M (2)	NWB file
Genotype (files)	Sst-IRES-Cre/wt;Ai32(RCL-ChR2(H134R)_EYFP)/wt (2)	NWB file
Strain (files)		NWB file
Age (files)	P144D (1), P163D (1)	NWB file
Lab / institution	/ Allen Institute for Neural Dynamics (2)	NWB file
Units (total)	7,490	NWB file
Units per file	median 3745 (range 3558–3932)	NWB file
Spikes (sum of spike_times_index lengths)		NWB file
How region was assigned to units	electrodes.location at the electrode (among the unit's units.electrodes) nearest to units.depth/estimated_x (file has no peak-channel column)	Computed by this export
Distinct region labels	79	NWB file
units.default_qc values	True: 4,542, False: 2,948	NWB file
units.decoder_label values	sua: 2,644, mua: 2,601, noise: 2,245	NWB file

The meaning of each units column, including any quality label, is quoted from the NWB file in section 7.



3. NWB files

File (session_id)	Subject	Sex	Age	Size (GB)	Units	Electrodes	Electrode groups	Interval tables (rows)
sub-767925_ses-ecephys-767925-2025-02-13-11-40-19_ecephys	767925	M	P144D	17.66	3932	2880	ProbeA, ProbeB, ProbeC, ProbeD, ProbeE, ProbeF	40 hz pulse train_presentations(149); 5 hz pulse train_presentations(149); behavior_presentations(3598); post_receptive_field_mapping_presentations(1944); post_surround_suppression_mapping_presentations(1500); pre_receptive_field_mapping_presentations(1944); pre_surround_suppression_mapping_presentations(1500); raised_cosine_presentations(150); spontaneous_presentations(6599)
sub-773206_ses-ecephys-773206-2025-04-01-10-47-58_ecephys	773206	M	P163D	16.80	3558	2880	ProbeA, ProbeB, ProbeC, ProbeD, ProbeE, ProbeF	40 hz pulse train_presentations(150); 5 hz pulse train_presentations(149); behavior_presentations(3597); post_receptive_field_mapping_presentations(1944); post_surround_suppression_mapping_presentations(1500); pre_receptive_field_mapping_presentations(1944); pre_surround_suppression_mapping_presentations(1500); raised_cosine_presentations(150); spontaneous_presentations(6598)

Other data in each file (not exported)

File (session_id)	Acquisition	Processing modules: interfaces
sub-767925_ses-ecephys-767925-2025-02-13-11-40-19_ecephys	autorewarded, licks, raw_running_wheel_rotation, reward_volume, running_wheel_signal_voltage, running_wheel_supply_voltage	ecephys:LFP; eye_tracking:likely_blink_times/corneal_reflection/ellipse/pupil; running:running_speed/running_wheel_rotation; stimulus:timestamps
sub-773206_ses-ecephys-773206-2025-04-01-10-47-58_ecephys	autorewarded, licks, raw_running_wheel_rotation, reward_volume, running_wheel_signal_voltage, running_wheel_supply_voltage	ecephys:LFP; eye_tracking:likely_blink_times/corneal_reflection/ellipse/pupil; running:running_speed/running_wheel_rotation; stimulus:timestamps

4. Units per region (Computed by this export)

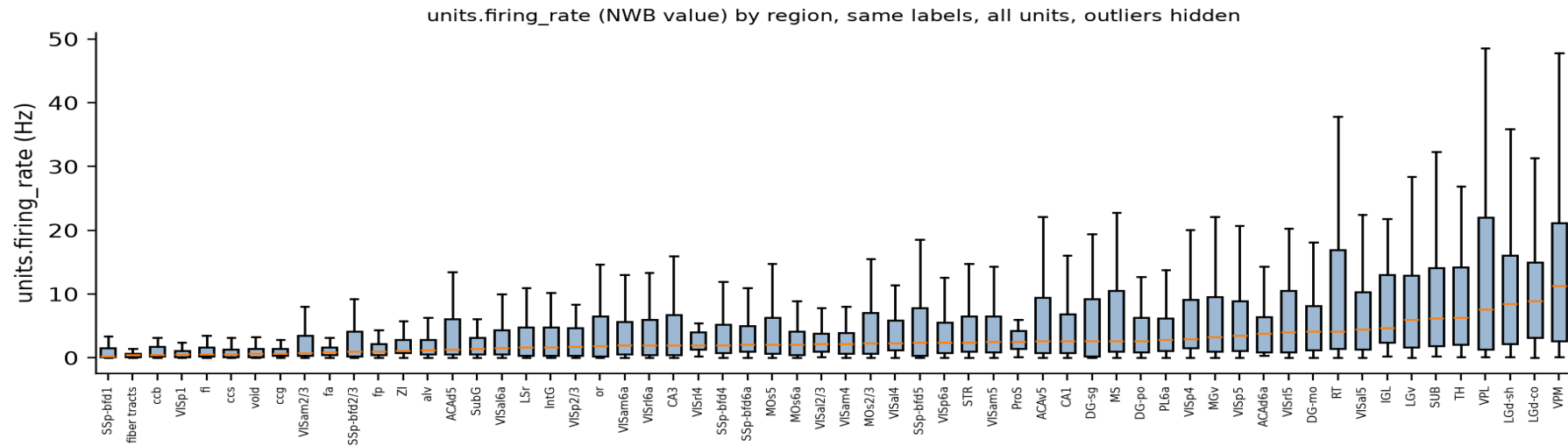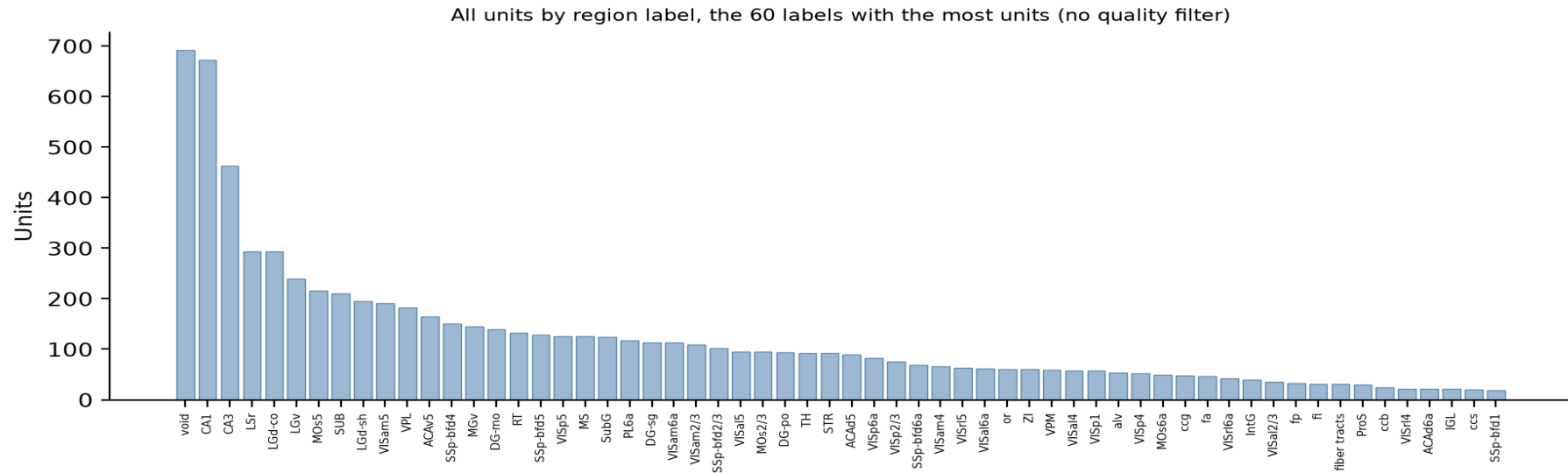

5. All region labels (Computed by this export)

One row per region label exactly as stored in the NWB files, all units (no quality filter), sorted by unit count. Firing-rate columns summarise the NWB units.firing_rate values where the dataset provides them.

Region label	Units	Files	Subjects	firing_rate mean (Hz)	firing_rate median (Hz)
void	692	2	2	1.18	0.63
CA1	672	2	2	5.67	2.59
CA3	463	2	2	5.42	1.98
LSr	293	2	2	3.69	1.66
LGd-co	293	2	2	10.69	8.91
LGv	239	2	2	8.90	5.97
MOs5	216	2	2	4.88	2.09
SUB	211	2	2	9.10	6.21
LGd-sh	195	2	2	11.41	8.43
VISam5	191	2	2	4.77	2.49
VPL	183	1	1	12.68	7.57
ACAv5	164	1	1	6.34	2.57
SSp-bfd4	151	2	2	4.40	2.00
MGv	145	1	1	6.84	3.23
DG-mo	140	2	2	5.46	4.11
RT	132	2	2	11.55	4.18
SSp-bfd5	128	2	2	5.41	2.40
VIsp5	125	2	2	6.35	3.45
MS	125	1	1	7.89	2.67
SubG	124	1	1	3.68	1.38
PL6a	118	1	1	4.19	2.83
DG-sg	113	2	2	7.25	2.64
VISam6a	113	2	2	4.06	1.93
VISam2/3	109	2	2	2.88	0.82
SSp-bfd2/3	102	2	2	3.29	0.97
VISal5	95	2	2	7.15	4.50
MOs2/3	95	2	2	4.96	2.27
DG-po	94	2	2	4.85	2.67
TH	92	2	2	9.63	6.34
STR	92	2	2	4.80	2.44
ACAd5	90	1	1	4.33	1.34
VIsp6a	82	2	2	4.26	2.42
VIsp2/3	76	2	2	3.20	1.75
SSp-bfd6a	69	2	2	4.10	2.03
VISam4	66	2	2	3.57	2.24
VISrl5	63	1	1	7.77	4.02
VISal6a	62	2	2	4.07	1.58
or	60	2	2	4.73	1.90
ZI	60	1	1	3.55	1.06
VPM	59	1	1	15.57	11.25

Region label	Units	Files	Subjects	firing_rate mean (Hz)	firing_rate median (Hz)
VISal4	58	2	2	4.74	2.28
VIsp1	58	2	2	1.23	0.51
alv	53	2	2	2.23	1.22
VIsp4	52	2	2	6.84	2.98
MOs6a	50	1	1	4.14	2.12
ccg	48	2	2	1.27	0.71
fa	46	1	1	1.67	0.92
VISrl6a	42	1	1	4.15	1.96
IntG	40	2	2	4.84	1.68
VISal2/3	35	2	2	3.62	2.15
fp	33	2	2	1.65	1.03
fi	31	1	1	1.15	0.53
fiber tracts	31	2	2	1.39	0.50
ProS	30	2	2	3.83	2.51
ccb	25	1	1	1.05	0.50
VISrl4	22	1	1	3.50	1.99
ACAAd6a	21	1	1	4.89	3.84
IGL	21	1	1	8.09	4.69
ccs	20	2	2	1.05	0.58
SSp-bfd1	19	2	2	0.89	0.27
VISrl2/3	17	1	1	5.87	4.47
VIsp6b	17	2	2	9.28	4.38
VISal6b	16	2	2	2.81	0.70
MGd	14	1	1	9.24	4.96
ACAv2/3	11	1	1	1.85	0.72
VISrl6b	10	1	1	2.99	1.74
VISam6b	10	1	1	3.23	2.09
CA2	9	2	2	3.97	3.36
em	8	1	1	4.20	3.69
SSp-bfd6b	8	2	2	5.70	2.49
VISam1	8	1	1	3.58	2.43
dhc	8	1	1	3.81	2.79
PL6b	7	1	1	2.66	2.24
bsc	7	1	1	1.96	0.73
VISal1	6	2	2	1.40	0.89
cing	4	1	1	0.51	0.47
HPF	1	1	1	2.74	2.74
MOs1	1	1	1	0.73	0.73
root	1	1	1	0.29	0.29

6. Intervals tables (stimuli, trials, epochs)

Each NWB intervals table, with the number of files containing it, total rows, and its columns as stored.

Table	Files	Rows	Columns
40_hz_pulse_train_presentations	2	299	start_time, stop_time, condition_num, level, condition, stim_name, duration
5_hz_pulse_train_presentations	2	298	start_time, stop_time, condition_num, level, condition, stim_name, duration
behavior_presentations	2	7,195	start_time, stop_time, duration, image_index, image_name, image_set, omitted, stim_name
post_receptive_field_mapping_presentations	2	3,888	start_time, stop_time, duration, stim_name, contrast, temporal_frequency, spatial_frequency, orientation, pos_x, pos_y
post_surround_suppression_mapping_presentations	2	3,000	start_time, stop_time, duration, stim_name, contrast, temporal_frequency, spatial_frequency, orientation, size, pos_x, pos_y
pre_receptive_field_mapping_presentations	2	3,888	start_time, stop_time, duration, stim_name, contrast, temporal_frequency, spatial_frequency, orientation, pos_x, pos_y
pre_surround_suppression_mapping_presentations	2	3,000	start_time, stop_time, duration, stim_name, contrast, temporal_frequency, spatial_frequency, orientation, size, pos_x, pos_y
raised_cosine_presentations	2	300	start_time, stop_time, condition_num, level, condition, stim_name, duration
spontaneous_presentations	2	13,197	start_time, stop_time, duration, stim_name

7. Column and signal descriptions from the NWB files

Quoted verbatim from sub-767925_ses-ecephys-767925-2025-02-13-11-40-19_ecephys.nwb; all files in the Dandiset share this structure.

Location in NWB file	Name	Description (as written in the NWB file)
electrodes	(table)	metadata about extracellular electrodes
electrodes	location	the location of channel within the subject e.g. brain region
electrodes	group	a reference to the ElectrodeGroup this electrode is a part of
electrodes	group_name	the name of the ElectrodeGroup this electrode is a part of
electrodes	channel_name	unique channel reference
electrodes	rel_y	rel_y
electrodes	rel_x	rel_x
electrodes	inter_sample_shift	no description
electrodes	x	ccf x coordinate
electrodes	y	ccf y coordinate
electrodes	z	ccf z coordinate
units	(table)	Units from Kilosort4
units	spike_times	the spike times for each unit
units	electrodes	the electrodes that each spike unit came from
units	waveform_mean	the spike waveform mean for each spike unit
units	waveform_sd	the spike waveform standard deviation for each spike unit
units	unit_name	Unique reference for each unit.
units	sync_spike_4	No description.
units	presence_ratio	No description.
units	amplitude_cv_range	No description.
units	isi_violations_ratio	No description.
units	nn_miss_rate	No description.
units	estimated_y	No description.
units	repolarization_slope	No description.
units	d_prime	No description.
units	num_positive_peaks	No description.
units	drift_mad	No description.
units	amplitude	No description.
units	estimated_z	No description.
units	velocity_above	No description.
units	nn_hit_rate	No description.
units	device_name	No description.
units	firing_range	No description.
units	amplitude_cv_median	No description.
units	exp_decay	No description.
units	estimated_x	No description.
units	rp_violations	No description.
units	velocity_below	No description.
units	amplitude_cutoff	No description.
units	recovery_slope	No description.

Location in NWB file	Name	Description (as written in the NWB file)
units	drift_std	No description.
units	snr	The signal-to-noise ratio of the unit.
units	isolation_distance	No description.
units	num_spikes	No description.
units	depth	No description.
units	sync_spike_2	No description.
units	ks_unit_id	No description.
units	default_qc	No description.
units	l_ratio	No description.
units	decoder_label	No description.
units	shank	No description.
units	drift_ptp	No description.
units	original_cluster_id	No description.
units	sync_spike_8	No description.
units	isi_violations_count	No description.
units	num_negative_peaks	No description.
units	decoder_probability	No description.
units	spread	No description.
units	half_width	No description.
units	firing_rate	Number of spikes per unit of time.
units	amplitude_median	No description.
units	sliding_rp_violation	No description.
units	silhouette	No description.
units	peak_trough_ratio	No description.
units	rp_contamination	No description.
units	peak_to_valley	The duration between the negative and the positive peaks computed on the maximum channel.
intervals/40 hz pulse train_presentations	(table)	Presentation times and stimuli details for '40 hz pulse train' stimuli. Note: image_name references control_description in stimulus/templates
intervals/40 hz pulse train_presentations	start_time	Start time of epoch, in seconds
intervals/40 hz pulse train_presentations	stop_time	Stop time of epoch, in seconds
intervals/40 hz pulse train_presentations	condition_num	Description of condition_num
intervals/40 hz pulse train_presentations	level	Description of level
intervals/40 hz pulse train_presentations	condition	Description of condition
intervals/40 hz pulse train_presentations	stim_name	Description of stim_name
intervals/40 hz pulse train_presentations	duration	Description of duration
intervals/40 hz pulse train_presentations	timeseries	index into a TimeSeries object
intervals/5 hz pulse train_presentations	(table)	Presentation times and stimuli details for '5 hz pulse train' stimuli. Note: image_name references control_description in stimulus/templates

Location in NWB file	Name	Description (as written in the NWB file)
intervals/5 hz pulse train_presentations	start_time	Start time of epoch, in seconds
intervals/5 hz pulse train_presentations	stop_time	Stop time of epoch, in seconds
intervals/5 hz pulse train_presentations	condition_num	Description of condition_num
intervals/5 hz pulse train_presentations	level	Description of level
intervals/5 hz pulse train_presentations	condition	Description of condition
intervals/5 hz pulse train_presentations	stim_name	Description of stim_name
intervals/5 hz pulse train_presentations	duration	Description of duration
intervals/5 hz pulse train_presentations	timeseries	index into a TimeSeries object
intervals/behavior_presentations	(table)	Presentation times and stimuli details for 'behavior' stimuli. Note: image_name references control_description in stimulus/templates
intervals/behavior_presentations	start_time	Start time of epoch, in seconds
intervals/behavior_presentations	stop_time	Stop time of epoch, in seconds
intervals/behavior_presentations	duration	Description of duration
intervals/behavior_presentations	image_name	Description of image_name
intervals/behavior_presentations	image_set	Description of image_set
intervals/behavior_presentations	omitted	Description of omitted
intervals/behavior_presentations	stim_name	Description of stim_name
intervals/behavior_presentations	timeseries	index into a TimeSeries object
intervals/injection_times	(table)	Intervals from ../data/interval_metadata/injection_times.csv
intervals/injection_times	start_time	Start time of epoch, in seconds
intervals/injection_times	stop_time	Stop time of epoch, in seconds
intervals/injection_times	lims_session_id	lims_session_id column from CSV
intervals/injection_times	session_id	session_id column from CSV
intervals/injection_times	injection_type	injection_type column from CSV
intervals/post_receptive_field_mapping_presentations	(table)	Presentation times and stimuli details for 'post_receptive_field_mapping' stimuli. Note: image_name references control_description in stimulus/templates
intervals/post_receptive_field_mapping_presentations	start_time	Start time of epoch, in seconds
intervals/post_receptive_field_mapping_presentations	stop_time	Stop time of epoch, in seconds
intervals/post_receptive_field_mapping_presentations	duration	Description of duration
intervals/post_receptive_field_mapping_presentations	stim_name	Description of stim_name
intervals/post_receptive_field_mapping_presentations	contrast	Description of contrast
intervals/post_receptive_field_mapping_presentations	temporal_frequency	Description of temporal_frequency
intervals/post_receptive_field_mapping_presentations	spatial_frequency	Description of spatial_frequency
intervals/post_receptive_field_mapping_presentations	orientation	Description of orientation
intervals/post_receptive_field_mapping_presentations	pos_x	Description of pos_x
intervals/post_receptive_field_mapping_presentations	pos_y	Description of pos_y

Location in NWB file	Name	Description (as written in the NWB file)
intervals/post_receptive_field_mapping_presentations	timeseries	index into a TimeSeries object
intervals/post_surround_suppression_mapping_presentations	(table)	Presentation times and stimuli details for 'post_surround_suppression_mapping' stimuli. Note: image_name references control_description in stimulus/templates
intervals/post_surround_suppression_mapping_presentations	start_time	Start time of epoch, in seconds
intervals/post_surround_suppression_mapping_presentations	stop_time	Stop time of epoch, in seconds
intervals/post_surround_suppression_mapping_presentations	duration	Description of duration
intervals/post_surround_suppression_mapping_presentations	stim_name	Description of stim_name
intervals/post_surround_suppression_mapping_presentations	contrast	Description of contrast
intervals/post_surround_suppression_mapping_presentations	temporal_frequency	Description of temporal_frequency
intervals/post_surround_suppression_mapping_presentations	spatial_frequency	Description of spatial_frequency
intervals/post_surround_suppression_mapping_presentations	orientation	Description of orientation
intervals/post_surround_suppression_mapping_presentations	size	Description of size
intervals/post_surround_suppression_mapping_presentations	pos_x	Description of pos_x
intervals/post_surround_suppression_mapping_presentations	pos_y	Description of pos_y
intervals/post_surround_suppression_mapping_presentations	timeseries	index into a TimeSeries object
intervals/pre_receptive_field_mapping_presentations	(table)	Presentation times and stimuli details for 'pre_receptive_field_mapping' stimuli. Note: image_name references control_description in stimulus/templates
intervals/pre_receptive_field_mapping_presentations	start_time	Start time of epoch, in seconds
intervals/pre_receptive_field_mapping_presentations	stop_time	Stop time of epoch, in seconds
intervals/pre_receptive_field_mapping_presentations	duration	Description of duration
intervals/pre_receptive_field_mapping_presentations	stim_name	Description of stim_name
intervals/pre_receptive_field_mapping_presentations	contrast	Description of contrast
intervals/pre_receptive_field_mapping_presentations	temporal_frequency	Description of temporal_frequency
intervals/pre_receptive_field_mapping_presentations	spatial_frequency	Description of spatial_frequency
intervals/pre_receptive_field_mapping_presentations	orientation	Description of orientation
intervals/pre_receptive_field_mapping_presentations	pos_x	Description of pos_x

Location in NWB file	Name	Description (as written in the NWB file)
intervals/pre_receptive_field_mapping_presentations	pos_y	Description of pos_y
intervals/pre_receptive_field_mapping_presentations	timeseries	index into a TimeSeries object
intervals/pre_surround_suppression_mapping_presentations	(table)	Presentation times and stimuli details for 'pre_surround_suppression_mapping' stimuli. Note: image_name references control_description in stimulus/templates
intervals/pre_surround_suppression_mapping_presentations	start_time	Start time of epoch, in seconds
intervals/pre_surround_suppression_mapping_presentations	stop_time	Stop time of epoch, in seconds
intervals/pre_surround_suppression_mapping_presentations	duration	Description of duration
intervals/pre_surround_suppression_mapping_presentations	stim_name	Description of stim_name
intervals/pre_surround_suppression_mapping_presentations	contrast	Description of contrast
intervals/pre_surround_suppression_mapping_presentations	temporal_frequency	Description of temporal_frequency
intervals/pre_surround_suppression_mapping_presentations	spatial_frequency	Description of spatial_frequency
intervals/pre_surround_suppression_mapping_presentations	orientation	Description of orientation
intervals/pre_surround_suppression_mapping_presentations	size	Description of size
intervals/pre_surround_suppression_mapping_presentations	pos_x	Description of pos_x
intervals/pre_surround_suppression_mapping_presentations	pos_y	Description of pos_y
intervals/pre_surround_suppression_mapping_presentations	timeseries	index into a TimeSeries object
intervals/raised_cosine_presentations	(table)	Presentation times and stimuli details for 'raised_cosine' stimuli. Note: image_name references control_description in stimulus/templates
intervals/raised_cosine_presentations	start_time	Start time of epoch, in seconds
intervals/raised_cosine_presentations	stop_time	Stop time of epoch, in seconds
intervals/raised_cosine_presentations	condition_num	Description of condition_num
intervals/raised_cosine_presentations	level	Description of level
intervals/raised_cosine_presentations	condition	Description of condition
intervals/raised_cosine_presentations	stim_name	Description of stim_name
intervals/raised_cosine_presentations	duration	Description of duration
intervals/raised_cosine_presentations	timeseries	index into a TimeSeries object
intervals/spontaneous_presentations	(table)	Presentation times and stimuli details for 'spontaneous' stimuli. Note: image_name references control_description in stimulus/templates
intervals/spontaneous_presentations	start_time	Start time of epoch, in seconds
intervals/spontaneous_presentations	stop_time	Stop time of epoch, in seconds
intervals/spontaneous_presentations	duration	Description of duration
intervals/spontaneous_presentations	stim_name	Description of stim_name
intervals/spontaneous_presentations	timeseries	index into a TimeSeries object
acquisition	autorewarded	Boolean flag indicating whether reward was auto-delivered
acquisition	licks	Lick events

Location in NWB file	Name	Description (as written in the NWB file)
acquisition	raw_running_wheel_rotation	no description
acquisition	reward_volume	Volume of reward delivered
acquisition	running_wheel_signal_voltage	no description
acquisition	running_wheel_supply_voltage	no description
processing/ecephys	(module)	Intermediate data from extracellular electrophysiology recordings, e.g., LFP.
processing/ecephys/LFP	ElectricalSeriesProbeA-LFP	LFP voltage traces from ProbeA
processing/ecephys/LFP	ElectricalSeriesProbeB-LFP	LFP voltage traces from ProbeB
processing/ecephys/LFP	ElectricalSeriesProbeC-LFP	LFP voltage traces from ProbeC
processing/ecephys/LFP	ElectricalSeriesProbeD-LFP	LFP voltage traces from ProbeD
processing/ecephys/LFP	ElectricalSeriesProbeE-LFP	LFP voltage traces from ProbeE
processing/ecephys/LFP	ElectricalSeriesProbeF-LFP	LFP voltage traces from ProbeF
processing/eye_tracking	(module)	eye tracking data
processing/eye_tracking	likely_blink_times	no description
processing/eye_tracking/corneal_reflection	(table)	
processing/eye_tracking/corneal_reflection	reference_frame	no description
processing/eye_tracking/corneal_reflection	data_x	no description
processing/eye_tracking/corneal_reflection	data_y	no description
processing/eye_tracking/corneal_reflection	area	no description
processing/eye_tracking/corneal_reflection	area_raw	no description
processing/eye_tracking/corneal_reflection	width	no description
processing/eye_tracking/corneal_reflection	height	no description
processing/eye_tracking/corneal_reflection	angle	no description
processing/eye_tracking/corneal_reflection	timestamps	no description
processing/eye_tracking/ellipse	(table)	
processing/eye_tracking/ellipse	reference_frame	no description
processing/eye_tracking/ellipse	data_x	no description
processing/eye_tracking/ellipse	data_y	no description
processing/eye_tracking/ellipse	area	no description
processing/eye_tracking/ellipse	area_raw	no description
processing/eye_tracking/ellipse	width	no description
processing/eye_tracking/ellipse	height	no description
processing/eye_tracking/ellipse	angle	no description
processing/eye_tracking/ellipse	timestamps	no description
processing/eye_tracking/pupil	(table)	
processing/eye_tracking/pupil	reference_frame	no description

Location in NWB file	Name	Description (as written in the NWB file)
processing/eye_tracking/pupil	data_x	no description
processing/eye_tracking/pupil	data_y	no description
processing/eye_tracking/pupil	area	no description
processing/eye_tracking/pupil	area_raw	no description
processing/eye_tracking/pupil	width	no description
processing/eye_tracking/pupil	height	no description
processing/eye_tracking/pupil	angle	no description
processing/eye_tracking/pupil	timestamps	no description
processing/running	(module)	running speed data
processing/running	running_speed	no description
processing/running	running_wheel_rotation	no description
processing/stimulus	(module)	Stimulus Times processing
processing/stimulus	timestamps	no description

8. Exported files

File	Rows	Columns
tables/sessions.csv	2	session_id, asset_path, size_gb, subject_id, species, sex, age, genotype, strain, session_start_time, session_description, institution, lab, n_units, n_electrodes, electrode_groups, region_source, regions, intervals, acquisition, processing
tables/units.parquet	7,490	session_id, unit_id, unit_row, region, unit_name, sync_spike_4, presence_ratio, amplitude_cv_range, isi_violations_ratio, nn_miss_rate, estimated_y, repolarization_slope, d_prime, num_positive_peaks, drift_mad, amplitude, estimated_z, velocity_above, nn_hit_rate, device_name, firing_range, amplitude_cv_median, exp_decay, estimated_x, rp_violations, velocity_below, amplitude_cutoff, recovery_slope, drift_std, snr, isolation_distance, num_spikes, depth, sync_spike_2, ks_unit_id, default_qc, l_ratio, decoder_label, shank, drift_ptp, original_cluster_id, sync_spike_8, isi_violations_count, num_negative_peaks, decoder_probability, spread, half_width, firing_rate, amplitude_median, sliding_rp_violation, silhouette, peak_trough_ratio, rp_contamination, peak_to_valley, subject_id
tables/region_summary.csv	79	region, n_units, n_sessions, n_subjects, firing_rate_mean_hz, firing_rate_median_hz
tables/intervals_summary.csv	9	table, n_sessions, n_rows, columns
intervals/<table>.parquet	one per intervals table	session_id + the table's columns
tables/dataset_metadata.json, nwb_file_metadata.json, nwb_descriptions.csv		published metadata quoted in sections 1 and 7